Double Flash Fusion

PROTOCOL DESCRIPTION DOCUMENT

IMOPA, Nancy

2026

PREPARING THE EXPERIMENTAL SET UP (BEFORE THE SUBJECT ARRIVES)

Note: The subject must be given an information form which informs them about their rights at least 24 hours before the start of the experiment.

- Prepare external electrodes (place ring stickers)
- Start computers: acq, stim, eyetracker
- Set subject ID in the 'sub00.csv' format on the stimulation scripts:
DFF_EXP/TempoSpacioTask_recording.mat
- Make sure the refresh rate is at 144hz in the monitor windows setting and in the script
- Start the eyetracker
- Pre-fill Fiche de renseignements
- Fill in a syringe with gel

PREPARING THE SUBJECT

- Ask the participant to sign the consent form
- Ask the subject if S/he needs to use the bathroom
- Let the subject be seated in the room and make sure s/he is comfortable.
- Measure the size of the head and choose the appropriate cap.
- Measure the position of Cz electrode and mark it
- Place the EEG cap so that Cz is at this point and such that electrodes are symmetric L/R & back/front.
- Place external electrodes:

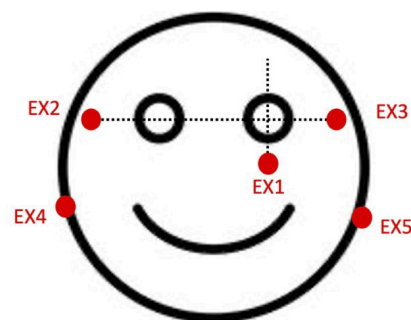
Fig 1. Place of external electrodes

- First run the [DFF_EXP/TempoSpacioTask_behav.mat](#) for then [DFF_EXP/TempoSpacioTask_recording.mat](#)

COLLECTING THE DATA

Configuring the eye-tracker

- Turn on the eye link by plugging it in
- Turn on the eyeblink host computer.
- Make sure the link is closed.
- Make sure the selected set options are as below.
- Ensure the subject's distance to the display screen is 80 cm.
- Put the marker on the subject's forehead.
- Start the camera setup and modify the parameters until you see the pupil size ok sign is green.
- Modify the lens to see that the image is no longer blurry.
- Choose both eyes.
- Choose the auto threshold.
- The illuminator power must be 50 per cent.
- Turn off the lights to start experiment

EEG

- Start Biosemi ActiView and make sure the right config file is loaded (e.g. double check the sampling rate on the right panel)
- Check the impedance on Biosemi (tab « Electrode Offset ») and adjust some electrodes if needed (any > 25).
- Make sure the subject understands the task.
- In Actiview (acquisition computer) start a new recording. (keep the recording on pause)
- Naming convention: sub00.bdf

Finishing

- Feedback: Ask the subject about the experiment and note answers on the subject recording sheet.
- Offer the subjects to clean their hair
- Clean the electrodes
(Warning: do not let the electrodes touch any metallic parts)
- Thank the subject for taking part in the experiment.

Equipment

- Stimulation computer:
 - Windows: 10
 - MATLAB: 2023a
 - Psychtoolbox
 - Screen: BENQ XL2430-B

- EEG Acquisition computer :
 - o Software: ActiView902-win.vi
 - o Config: Default_DFF.cfg
- Eye tracker computer :
 - o Eyelink 1000 plus
 - o Sampling rate: 500 Hz
- BioSemi System 64 + 8 electrodes

PARADIGM DESCRIPTION AND TRIGGERS

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